



Half Screw, Half Wire...
 Introducing the **Scwire®** by **OrthoPediatrics**



The world of fracture care has entered a new age...


Every good and noble hero needs special weapons to do battle everyday. The unique engineering and design of the multipurpose Scwire® provides agility, compression and superior fixation.

The Scwire offers: percutaneous application, reverse cutting flutes (facilitating easy removal), simple instrumentation, cut-to-length capability and a compression nut providing the ability to compress fragments after drilling.

With its “dual thread” design and accompanying compression nut, the Scwire is designed to stabilize comminuted displaced small bone fractures of the upper and lower extremities.

The Scwire is a simple and effective solution for small fracture fixation.


2.5 – Mini, 10 and 20

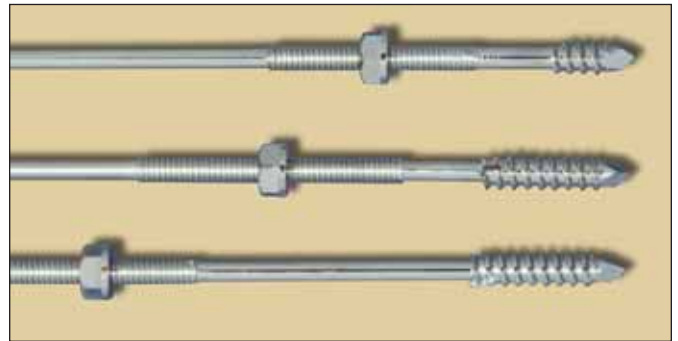

2.5 Anchor Threads



2.5 Compression Nut

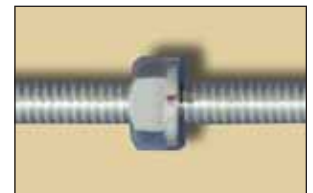

Indications

Percutaneous Compression Wire System for the fixation of osseous fragments of fractures of small bones, including:

Upper Extremity


- Finger Arthodesis
- Radial Head Fractures
- Bennett’s Fractures
- Olecranon Fractures
- Scaphoid Fractures

Lower Extremity

- Fore/Midfoot Fusions
- Metatarsal Fractures
- LisFranc Fractures
- Tillaux Fractures
- Ankle Fractures

Features and Benefits


- Simple instrumentation for a shorter learning curve
- Combined technology of a screw and wire reduces surgical steps and time
- Multiple sizes for a broader range for applications
- Compression capabilities provides a stable fracture reduction
- Small fracture/bone capability stabilizes fractures without losing reduction
- Small diameter allows for use near growth plates and for periarticular applications

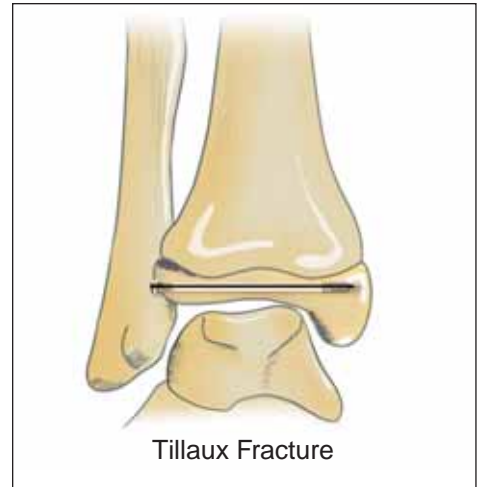

Scwire® Instrument Set

Instruments	Catalog Number
Socket Wrench	01-1070-003
Soft Tissue Protector	01-1070-004
Removal Tool	01-1070-005
Cannulated Cutter	01-1070-006
Counter Sink	01-1070-007
Case / Tray	01-1070-002

Complete Kit and Tray 01-1070-020

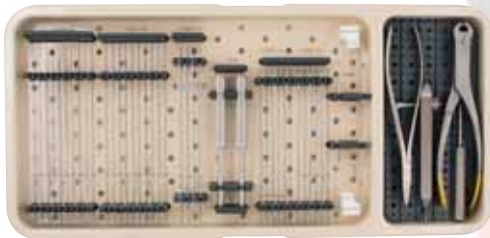

Scwire® 2.5 – Mini Implants

Stainless Steel

Catalog Number
00-1070-004



Scwire® 2.5 – 10 Implants

Stainless Steel

Catalog Number
00-1070-001



Additional Scwire® Components

Stainless Steel

Item	Catalog Number
2.5 FD (fusion device)	00-1070-003
2.5 Compression Nut	00-1070-005
2.5 Flat Washer	00-1070-006

Scwire® 2.5 – 20 Implants

Stainless Steel

Catalog Number
00-1070-002

